



All 3M Products Newsroom About 3M

Search



United States  [Change](#)

🏠 Infection Prevention

[By Profession](#)

[Products](#)

[CE Education](#)

[What's New](#)

[Contact Us](#)

[United States](#) > [Health Care](#) > [Infection Prevention](#) > [Breathability](#)

Not all N95s are alike.




Breathe easier with 3M.

3M™ surgical N95 respirators show lower pressure drop than all other competitive respirators tested.^{†*}

In general, a lower Delta-P translates to increased breathability~


Medical Respirator Differential Pressure (Delta-P) Study

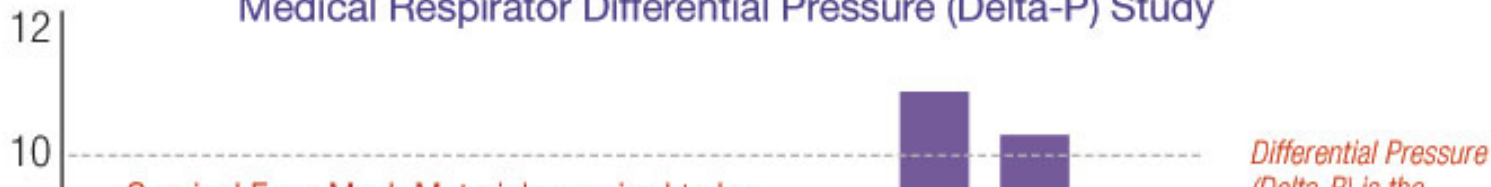

[Request Free Sample](#)



[1870](#) or [1870+](#)

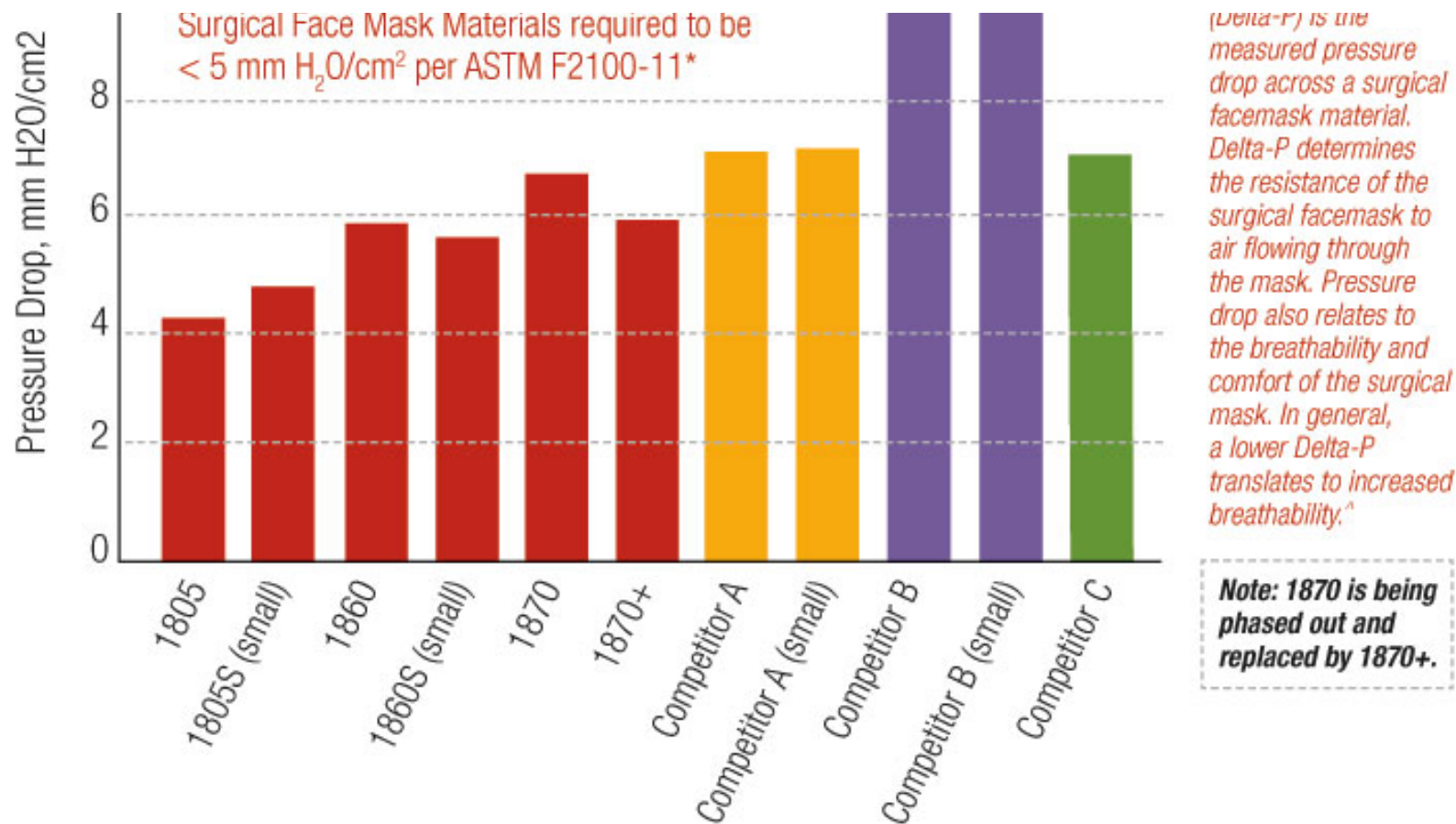

*Testing per MIL-Spec 36945C 4.4.1.2 using 4.9 cm² sample at 8 L/min. Source: Pressure drop performance in chart represents results from internal testing conducted on a sample lot of 3M products against a sample lot of commercially available product from three leading competitors in September 2012. Selection of respirator models, testing protocol, data generation and conclusions were reviewed and approved by an expert from the University of Minnesota. Statistically different per two sample T-tests (P-value <0.001). Results may vary.

*Performance requirement ASTM F2100-11 medical facemask materials classification.

[^]MIL-M 36945C 4.4.1.1.1 Method 1 Military Specifications: Surgical Mask, disposable (June 12, 1975) (Source – Guidance for Industry and FDA Staff: Surgical Mask-Premarket Notification [510 (k)] Submissions)



[1860](#) or [1860S](#)



[1805](#) or [1805S](#)



[Get Our Breathability Brochure](#)

Breathe easier knowing you have protection

Protect your patients, your health, your profession

In the health care setting, exposure to certain airborne particulates such as M. tuberculosis, those found in laser plumes and more may put professionals like you at risk. With advanced electrostatically-charged microfiber filter media, 3M™ Health Care respirators can help protect you while maintaining the comfort you have come to expect from the leader in respiratory protection technology.



Our unique manufacturing process injects a high level of electrostatic charge into microfibers arranged in an open formation that allows for the greater passage of air through the filter media.

3M Health Care Particulate Respirator and Surgical Masks are Engineered with proprietary Advanced Electrostatic Media – designed for ease of breathing. [See competitive study.](#)

Highly charged microfibers enhance the capture of airborne particles while helping you to breathe easier.


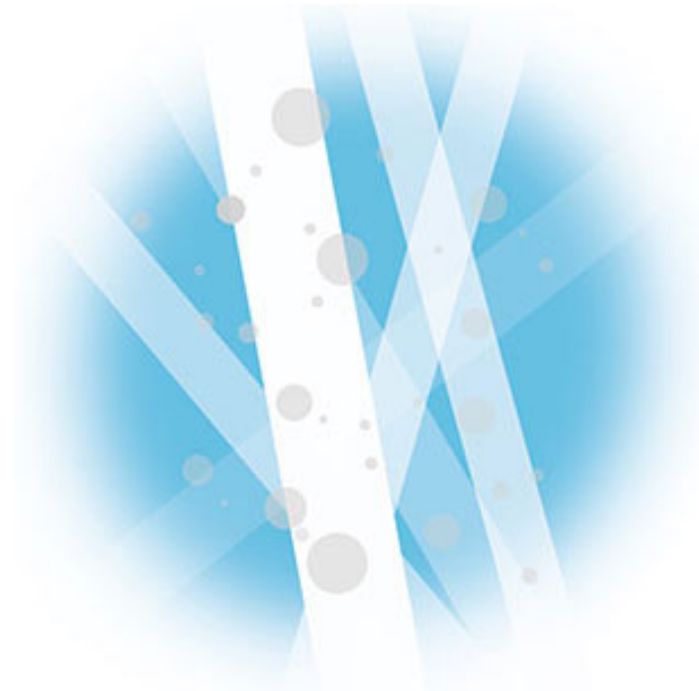


Illustration of 10,000x magnified electrostatically charged microfibers.

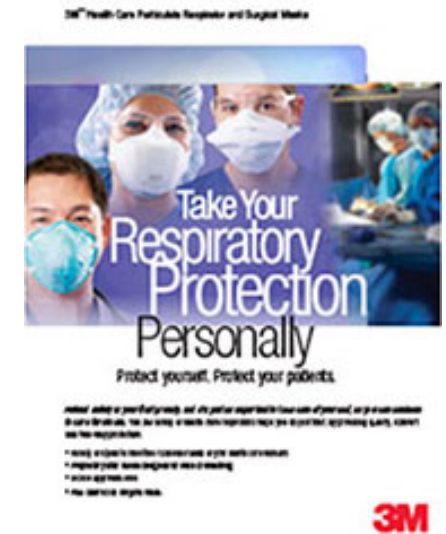

[Download 3M Health Care Respirator and Surgical Mask Family Brochure](#)

Infection Prevention: [About Us](#) | [Contact Us](#) | [Sign-Up for eMail](#)

3M: [News](#) | [Careers](#) | [Partners & Suppliers](#) | [Investor Relations](#) | [SDS Search](#) | [Transport Information Search](#) | [Help Center](#)

